



LightCOM

Lighthouse Consortium on
Obesity Management

Study name: **Lighthouse Consortium on Obesity Management**

Short name: **LightCOM**

SOP title: **Standard Operating Procedure for Continuous Glucose monitoring by Dexcom LightBAR**

SOP short title: **SOP on Dexcom**

SOP Version: **1**

Date: **22.10.2025**

Authored by:

Date:

Sabrina Wielsøe

Indhold

1. PURPOSE AND SCOPE.....	3
1.1 Set-up.....	3
1.2 Receiving and Data Upload	3


1. PURPOSE AND SCOPE

This standard operating procedure (SOP) defines the process for set-up and receiving Dexcom G7 monitors.

This SOP applies to all personnel working with Dexcom CGM in the LightBAR trial.

1.1 Set-up

1. Turn on the Receiver
2. Follow the instruction on the screen and pair with sensor
3. Once paired connect to PC via USB-cable
4. Log on to clarity.dexcom.eu
5. Choose Manage Receiver and blind

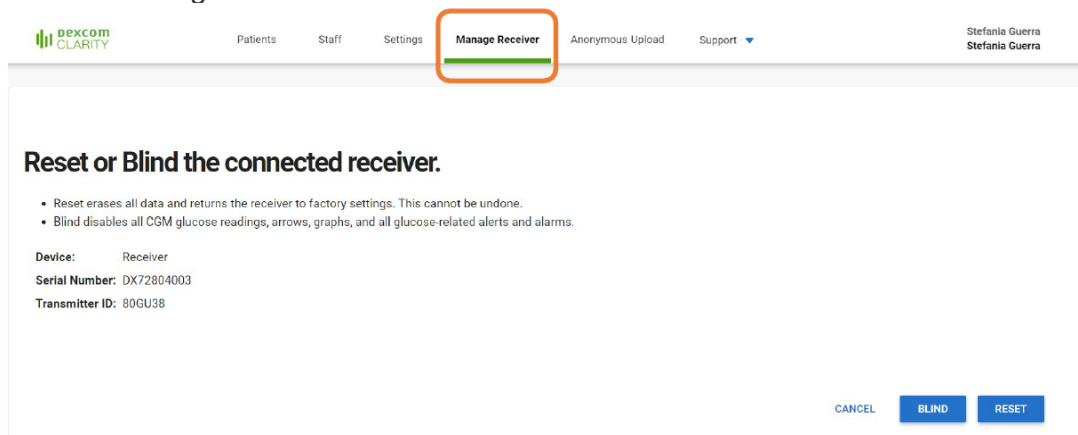

6. The Sensor can store data for 24 h. In order to transfer the data, the receiver and sensor should be within 6 m of each other.
We recommend the participant to have the receiver at the night table, the sensor and receiver will then connect during night at transfer data.

1.2 Receiving and Data Upload

7. Clean the receiver
8. Connect to PC via USB-cable
9. Log on to clarity.dexcom.eu
10. Choose Manage Receiver and unblind

11. Add Patient (or find Patient if follow-up)

Add new patient

Export all data

PATIENT NAME	DOB	PATIENT ID	LAST UPLOADED	DATA SHARING
guerra1, stefi	2 Oct 1983		20 May 2020	✔ On
Guerra 2, Stefi	2 Oct 1983	Study002		✖ Off
Pro, Stefi	2 Oct 2000		18 May 2020	✖ Off

Add new Patient

Add new Patient

First Name

Last Name

Date of Birth

DD

Month ▼

YYYY

Patient ID (optional) ?


Save

Cancel

12. Upload Data

Data upload


- To upload data, the receiver must be in **unblinded mode**. To upload data after a blinded session, please unblind the receiver through “manage receiver” as a first step.
- Click on Subject's “Name” in Patient List
- Click “Upload Data”
- To upload data, follow instructions on screen to connect receiver to computer with USB cable

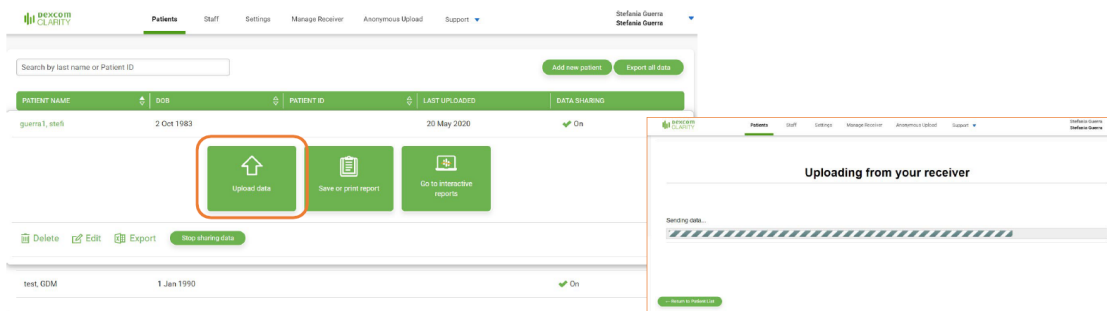

13. Data ensures

Ensure there is data for the desired period.

- “Go to interaction reports”
- Choose Daily pattern and check there is data for the desired dates.
- If anything is missing, please make a note in Open Clinica and Protocol Deviations

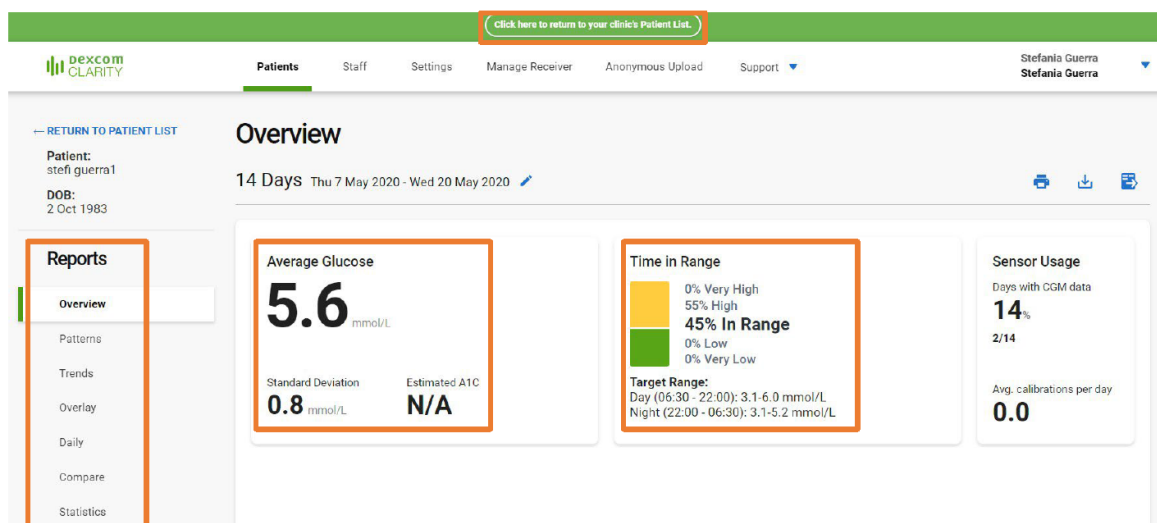

14. Reset and turn off Receiver:

- Go to Manage Receiver
- Reset
- Turn of the Receiver